

Syncerity BI UV-Vis

Scientific Deep-cooled Camera

OSD-SY-02

ELEMENTAL ANALYSIS
FLUORESCENCE
GRATINGS & OEM SPECTROMETERS
OPTICAL COMPONENTS
FORENSICS
PARTICLE CHARACTERIZATION
RAMAN
SPECTROSCOPIC ELLIPSOmetry
SPR IMAGING



UV-enhanced Sensitivity
with Ultra-high
Spectral Resolution!

Key Features and Benefits

- **2048 × 70 back-illuminated sensor**
Enable optimum spectral resolution
- **UV-Vis quantum efficiency enhancement**
60% QE at 250 nm, and 75% QE at 550 nm
- **Deep thermoelectric cooling**
−50°C for low dark current
- **Improved etaloning**
Ideal for Raman applications
- **16-bit digitization**
Provides wide dynamic range
- **Lifetime vacuum warranty**
Metal-sealed technology for permanent vacuum

Sensor Size 2048 × 70

Deep-cooled −50°C


Pixel Size 14 μm × 14 μm

Digitization 16 bit


Sample Applications

- Raman spectroscopy
- Microspectroscopy
- Plasma analysis
- UV-VIS-NIR photoluminescence
- Diffuse reflectance spectroscopy

Quantum Efficiency

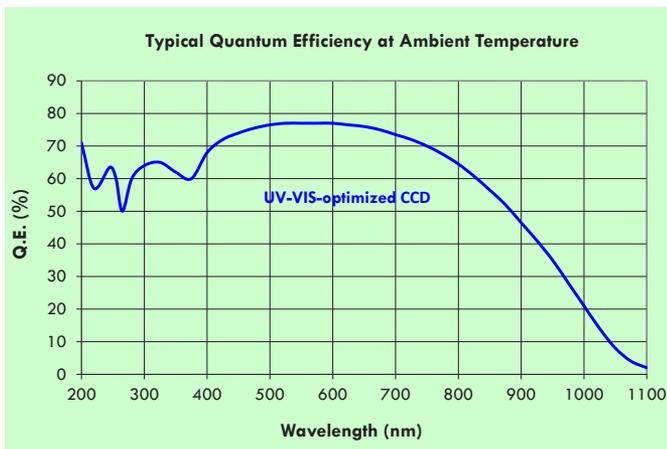


Suppressed Etaloning

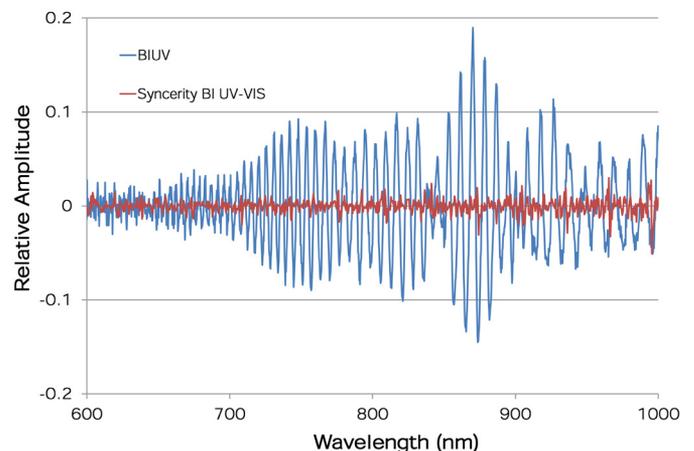

Syncerity™ BI UV-Vis Specifications

CCD Sensor Format	2048 × 70
Quantum efficiency at 20°C	63% at 245 nm; 64% at 300 nm; 68% at 400 nm 76% at 500 nm; 77% at 700 nm; 64% at 800 nm
Pixel size	14 μm × 14 μm
Image area	28.7 mm × 0.98 mm, 100% fill factor
Deep thermoelectric cooling	-50°C at +25°C ambient (-60°C at +25°C ambient on request) Yields low dark current suitable for most OEM and some research applications
Single pixel well capacity	50 000 e ⁻ /pixel (minimum); 60 000 e ⁻ /pixel (typical)
Serial register full well capacity	250 000 e ⁻ /pixel (minimum) 500 000 e ⁻ /pixel (typical output register saturation)
Scan rates	45 kHz and 500 kHz
Readout noise (at 45 kHz and at -50°C)¹	9 e ⁻ (typical) to 12 e ⁻ (maximum)
Readout noise (at 500 kHz and at -50°C)¹	20 e ⁻ (typical) to 25 e ⁻ (maximum)
Maximum spectral rate	20 Hz at 45 kHz scan rate 189 Hz at 500 kHz scan rate
Digitization	16-bit ADC
Dynamic range (typical for single pixel)²	55 500:1
Non-linearity (measured on each camera)	<0.15% (typical) at 45 kHz (0.4% maximum) <0.20% (typical) at 500 kHz (1% maximum)
Dark current at -50°C³ (Note that pixel size = 14 μm)	0.05 e ⁻ /pixel/s (typical)
Software-adjustable gains	2, 4, and 10 e ⁻ /count at -50°C
Environmental conditions	<ul style="list-style-type: none"> Operating temperature 0°C to 40°C ambient Relative humidity <70% (non-condensing) Storage temperature -25°C to 50°C
Weight	1.769 kg (3.90 lb)
Dimensions	See mechanical drawings
Power requirements AC/DC power supply (provided) Recommendation for OEM supplying camera to power directly:	90–264 VAC, 47–63 Hz <ul style="list-style-type: none"> Pin: +9 V, ± 5%, 6.44 A maximum Regulation: +8.55 V_{min}, +9 V_{typ}, +9.45 V_{max} Ripple & Noise: 200 mV_{pp} maximum
Minimum computer requirements	<ul style="list-style-type: none"> 3.0 GHz single core or 2.4 GHz multi-core processor 2 GB RAM 32-bit or 64-bit compatible 500 MB free hard disk space (additional disk space may be required depending on data-storage needs) USB 2.0 High-speed host controller capable of sustained rate of 40 MB/s Windows® (XP, Vista and 7)

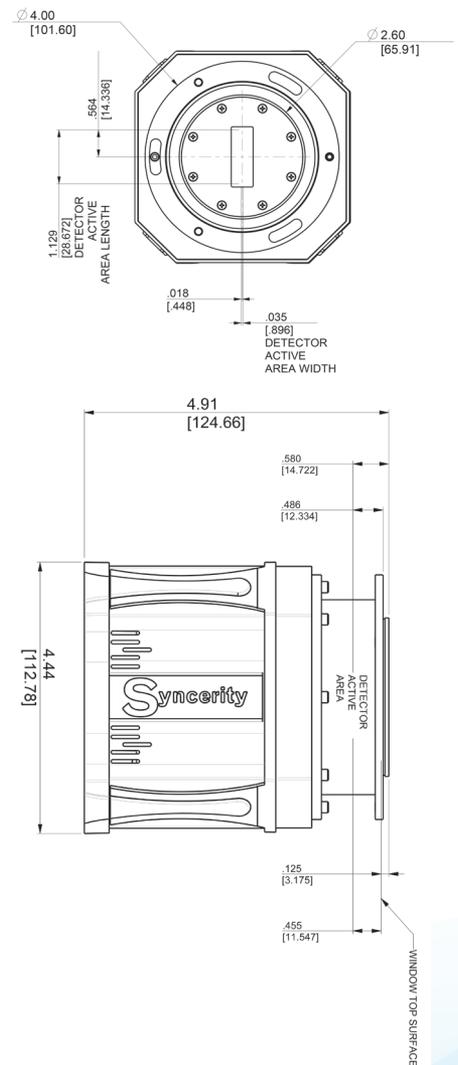

1. Entire system noise measured for a single pixel
2. Dynamic range is defined as Full Well/Readout Noise, measured at 45 kHz
3. Averaged over CCD area, but excluding any regions of blemishes.



Your partner for Spectroscopy Solutions



Spectrometers, Monochromators
Scientific Cameras, Software and
Custom Solutions



HORIBA
Scientific

info.sci@horiba.com

USA: +1 732 494 8660
UK: +44 (0)20 8204 8142
China: +86 (0)21 6289 6060

France: +33 (0)1 69 74 72 00
Italy: +39 2 5760 3050
Brazil: +55 (0)11 5545 1500

www.horiba.com/scientific

Germany: +49 (0)89 4623 17-0
Japan: +81 (0)3 6206 4721
Other: +1 732 494 8660

Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA